


Instruction Manual

HI 99163

Portable Waterproof Meat pH Meter

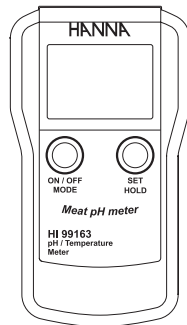

WARRANTY

All Hanna Instruments meters are warranted for two years against defects in workmanship and materials when used for their intended purpose and maintained according to instructions. The probes are warranted for a period of six months.

This warranty is limited to repair or replacement free of charge. Damage due to accidents, misuse, tampering or lack of prescribed maintenance are not covered.

If service is required, contact the dealer from whom you purchased the instrument. If under warranty, report the model number, date of purchase, serial number and the nature of the problem. If the repair is not covered by the warranty, you will be notified of the charges incurred. If the instrument is to be returned to Hanna Instruments, first obtain a Returned Goods Authorization Number from the Customer Service department and then send it with shipment costs prepaid. When shipping any instrument, make sure it is properly packaged for complete protection.

All rights are reserved. Reproduction in whole or in part is prohibited without the written consent of the copyright owner.

Dear Customer,
Thank you for choosing a Hanna Instruments Product.

Please read this instruction manual carefully before using the instrument.

If you need additional technical information, do not hesitate to e-mail us at tech@hannainst.com

This instrument is in compliance with the CE directives.

PRELIMINARY EXAMINATION

Remove the instrument from the packing material and examine it carefully to make sure that no damage has occurred during shipment. If noticeable damage is evident, notify your dealer.

Note: Save all packing material until you are sure that the instrument functions correctly. All defective items must be returned in the original packing together with the supplied accessories.

GENERAL DESCRIPTION

HI 99163 is a portable pH and temperature meter with a special probe, dedicated to the application in slaughter houses.


The FC 232D probe has been specially designed with a stainless steel blade tip for meat penetration.

Moreover, a built-in temperature sensor allows simultaneous pH and temperature readings and a mini amplifier renders measurements impervious to noise and electrical interferences.

The meter is supplied in a rugged carrying case, complete with:

- FC 232D amplified pH/temperature probe with FC 099 stainless steel blade, DIN connector and 1 m (3.3') cable
- HI 70004 pH 4.01 buffer, 20 mL sachet
- HI 70007 pH 7.01 buffer, 20 mL sachet
- HI 700630 Acid cleaning solution for meat grease and fats, 20 mL sachet (2 pcs.)
- batteries and instructions

LCD DESCRIPTION

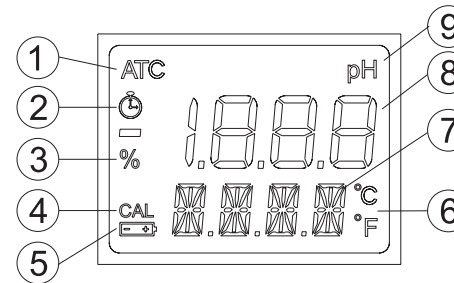

1. Automatic Temperature Compensation indicator
2. Stability indicator
3. Battery percentage
4. pH calibration indicator
5. Low battery indicator
6. Selectable temperature unit
7. Secondary display
8. Primary display

SPECIFICATIONS

Range	-2.00 to 16.00 pH -5.0 to 105.0°C / 23.0 to 221.0°F
Resolution	0.01 pH / 0.1°C / 0.1°F
Accuracy	±0.02 pH (@20°C/68°F) ±0.5°C up to 60°C; ±1°C outside ±1°F up to 140°F; ±2°F outside
Typical EMC Deviation	±0.02 pH ±0.2°C or ±0.4°F
Temp. Compensation	Automatic
pH Calibration	Automatic, 1 or 2 point with 2 sets of standardized buffers (4.01/7.01/10.01 or 4.01/6.86/9.18)
Probe (included)	FC 232D amplified pH/temperature probe
Battery Type	3 x 1.5V AA / IEC LR6
Battery Life	Approx. 1500 hours
Auto-off	After 8 minutes of non-use
Environment	0 to 50°C; RH 100%
Dimensions	150 x 80 x 36 mm (5.9 x 3.2 x 1.4")
Weight	210 g (7.4 oz.)

ACCESSORIES

- FC 232D Amplified pH probe with stainless steel blade tip, built-in temperature sensor, DIN connector and 1 m (3.3') cable
- FC 098 20 mm stainless steel blade
- FC 099 35 mm stainless steel blade
- HI 70004P pH 4.01 buffer solution, 20 mL sachet (25 pcs.)
- HI 7004L pH 4.01 buffer solution, 500 mL
- HI 70006P pH 6.86 buffer solution, 20 mL sachet (25 pcs.)
- HI 7006L pH 6.86 buffer solution, 500 mL
- HI 70007P pH 7.01 buffer solution, 20 mL sachet (25 pcs.)
- HI 7007L pH 7.01 buffer solution, 500 mL
- HI 70300M Storage solution, 230 mL bottle
- HI 70300L Storage solution, 500 mL bottle
- HI 7061M General purpose cleaning solution, 230 mL bottle
- HI 7061L General purpose cleaning solution, 500 mL bottle
- HI 700630P Acid cleaning solution for meat grease and fats, 20 mL sachet (25 pcs.)
- HI 70630M Acid cleaning solution for meat grease and fats, 230 mL bottle
- HI 70630L Acid cleaning solution for meat grease and fats, 500 mL bottle
- HI 70631M Alkaline cleaning sol. for meat grease and fats, 230 mL bottle
- HI 70631L Alkaline cleaning sol. for meat grease and fats, 500 mL bottle
- HI 70632M Cleaning & disinfection solution for blood products, 230 mL
- HI 70632L Cleaning & disinfection solution for blood products, 500 mL
- HI 710007 Blue shockproof rubber boot
- HI 710008 Orange shockproof rubber boot
- HI 721312 Hard carrying case

Hanna Instruments reserves the right to modify the design, construction and appearance of its products without advance notice.

OPERATIONAL GUIDE

To connect the probe

With the meter turned off, connect the **FC 232D** probe to the DIN socket on the top of the meter by aligning the pins and pushing in the plug. Tighten the nut to ensure a good connection. Remove the protective cap from the probe before taking any measurements.

To turn the meter ON and check the battery status

Press the ON/OFF/MODE button until the display lights up. At start-up, all the LCD segments are displayed for 1 second, then the percent indication of the remaining battery life is shown for another second (E.g. % 100 BATT). The meter then enters the normal measuring mode.

Note: If the display needs to be checked, keep the ON button pressed while turning the meter on. The meter will display all segments as long as the button is pressed.

To freeze the display

While in measurement mode, press the SET/HOLD button, HOLD appears on the secondary display and the reading will be frozen on the LCD (E.g. pH 5.73 HOLD). Press any button to return to normal mode.

To turn the meter OFF

While in normal measurement mode, press the ON/OFF/MODE button. OFF will appear on the secondary display. Release the button.

Note: The meter is provided with an acoustic signal feature, which can be disabled using the switch located in the battery compartment.

Note: When the meter detects the absence of a temperature probe at its input, the Automatic Temperature Compensation is turned off, and the meter uses a default value of 25°C (77°F) for the temperature measurement and compensation. In this condition, the secondary LCD shows 25.0°C (77.0°F) blinking. When a probe is connected, the meter automatically returns to the ATC mode, the ATC tag is turned on, and the temperature is shown on the secondary display.

Note: To clean the meter, do not use aggressive detergents. It is recommended to use water.

pH MEASUREMENT & CALIBRATION

- Make sure the meter has been calibrated before use.
- If the probe is dry, soak it in **HI 70300** storage solution for one hour to reactivate it.
- Insert the tip of the probe into the sample to be tested and wait until the □ symbol on the top left of the LCD is turned off.
- The LCD shows the pH value (automatically compensated for temperature) on the primary LCD, while the secondary LCD shows the temperature of the solution.
- After each measurement, it is recommended to perform a complete cleaning of the probe: use the suitable **Cleaning solutions** (see "Accessories" section) and clean the probe tip thoroughly to eliminate dirtiness and avoid cross-contamination.
- When not in use, switch the meter OFF and clean the probe with the **specific Cleaning solutions** (see "Accessories" section) before storing it.

pH calibration

For better accuracy, frequent calibration of the instrument is recommended. In addition, the instrument must be recalibrated whenever:

- a) The pH electrode is replaced.
 - b) After testing aggressive chemicals.
 - c) Where high accuracy is required.
 - d) At least once a month.
- From normal mode, press and hold the ON/OFF/MODE button until OFF on the secondary display is replaced by CAL. Release the button.
 - The LCD enters the calibration mode, displaying "pH 7.01 USE" (or "pH 6.86 USE" if the NIST buffer set was selected). After 1 second the meter activates the automatic buffer recognition feature. If a valid buffer is detected, then its value is shown on the primary display, and REC appears on the secondary LCD. If no valid buffer is detected, the meter keeps the USE indication active for 12 seconds, and then replaces it with WRNG, indicating that the sample being measured is not a valid buffer.

- For a single-point calibration with buffers pH 4.01, 9.18 or 10.01, the meter automatically accepts the calibration when the reading is stable; the meter will show on the primary display the accepted buffer, with the message "OK 1" on the secondary display, and an audible signal is produced. After 1 second the meter automatically returns to the normal measuring mode.
If a single-point calibration with buffers pH 7.01 or 6.86 is desired, then after the calibration point has been accepted the ON/OFF/MODE button must be pressed in order to return to the normal measuring mode. After the button is pressed, the meter shows "7.01" (or "6.86") - "OK 1", and an audible signal is produced. After 1 second, the meter automatically returns to the normal measuring mode.

Note: It is always recommended to carry out a two-point calibration for better accuracy.

- For a two-point calibration, place the probe in pH 7.01 (or pH 6.86) buffer. After the calibration point has been accepted, the "pH 4.01 USE" message appears. The message is held for 12 seconds, unless a valid buffer is recognized. If no valid buffer is recognized, then the WRNG message is shown. If a valid buffer (pH 4.01, pH 10.01 or pH 9.18) is detected, then the meter completes the calibration procedure. When the buffer is accepted, the LCD shows the accepted value with the "OK 2" message on the secondary display. The meter then returns to the normal measuring mode.

Note: When the calibration procedure is completed, the CAL tag is turned on.

To quit calibration and to reset to the default values

- After entering the calibration mode and before the first point is accepted, it is possible to quit the procedure and return to the last calibration data by pressing the ON/OFF/MODE button. The secondary LCD displays ESC for 1 second and the meter returns to normal mode.
- To reset to the default values and clear a previous calibration, press the SET/HOLD button after entering the calibration mode and before the first point is accepted.

The secondary LCD displays CLR for 1 second, the meter resets to the default calibration and the "CAL" tag on the LCD disappears.

METER SETUP

Setup mode allows the selection of the temperature unit and of the pH buffer set.

To enter the Setup mode, press and hold the ON/OFF/MODE button until CAL on the secondary display is replaced by TEMP and the current temperature unit (E.g. TEMP °C). Then:

- for °C/°F selection, use the SET/HOLD button. After the temperature unit has been selected, press ON/OFF/MODE to enter the buffer set selection mode; press ON/OFF/MODE twice to return to the normal measuring mode.
- to change the calibration buffer set, after setting the temperature unit, the meter will show the current buffer set: "pH 7.01 BUFF" (for 4.01/7.01/10.01) or "pH 6.86 BUFF" (for 4.01/6.86/9.18). Change the set with the SET/HOLD button, then press ON/OFF/MODE to return to normal mode.

BATTERY REPLACEMENT

The meter displays the remaining battery percentage when turned on. When the level is below 5%, the □ symbol on the bottom left of the LCD blinks to indicate a low battery condition. If the battery level is low enough to cause erroneous readings, the Battery Error Prevention System (BEPS) turns the meter off.

Unscrew the 4 screws located on the back of the meter and carefully replace the 3 AA batteries located in the battery compartment, while paying attention to their polarity. Reattach the back making sure that the gasket is in place and tighten the 4 screws to ensure a watertight seal.